

# Using Cholangioscopy to Investigate Anatomical Variation in Bile Duct Anatomy and Facilitate Stone Extraction through Selective Cannulation



**Dhavan Parikh, M.D.**

The Everett Clinic  
Everett, Washington, U.S.A.





## technique spotlight

### Patient History


A 76-year-old man was referred by an outside hospital and presented to our hospital with biliary colic. Computed tomography imaging at outside facility noted choledocholithiasis and cholelithiasis. Upon admission, the patient was noted to have elevated bilirubin (7.8 mg/dL) and alkaline phosphatase (346 U/L) levels. An initial ERCP with sphincterotomy was performed with clearance of the common bile duct. The patient underwent laparoscopic cholecystectomy the following day with intraoperative cholangiogram showing multiple filling defects in the common bile duct, presumably from dropped stones during the cholecystectomy.

### Procedure

Cannulation of the papilla was performed through the existing sphincterotomy and a 0.035" x 450cm straight Hydra Jagwire™ High Performance Guidewire was placed into the biliary tree. A cholangiogram showed multiple filling defects in the mid common bile duct (Figure 1). Sphincteroplasty up to 10mm was performed, followed by balloon sweeps that recovered numerous cholesterol stones. An occlusion cholangiogram was then performed to assess for clearance of the bile duct. This showed two filling defects adjacent to the distal common bile duct (Figure 2). A balloon sweep was then repeated but failed to clear the filling defects. The SpyGlass™ DS System SpyScope™ was passed over the existing guidewire to investigate further. Cholangioscopy was performed and uncovered atypical biliary anatomy with low medial insertion of the cystic duct (Figure 3). The guidewire was removed and the SpyScope was advanced into the cystic duct where two yellow stones were detected (Figure 4). The SpyScope

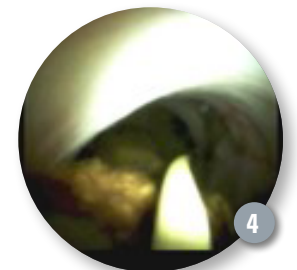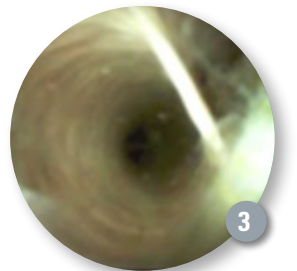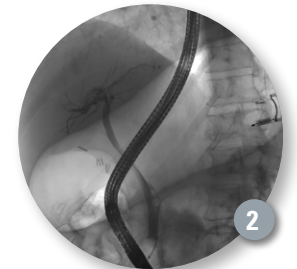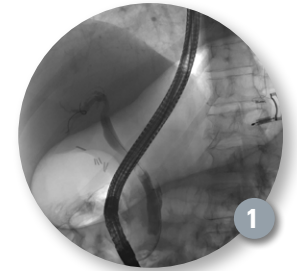

# Using Cholangioscopy to Investigate Anatomical Variation in Bile Duct Anatomy and Facilitate Stone Extraction through Selective Cannulation

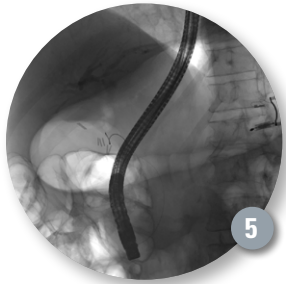


## technique spotlight

was then used to selectively cannulate the cystic duct, which was confirmed on fluoroscopy (Figure 5). The SpyScope™ was then exchanged for an extraction balloon and balloon sweeps of the cystic duct led to removal of both stones (Figure 6). Repeat occlusion cholangiogram through the main duct did not reveal any additional filling defects (Figure 7). The patient tolerated the procedure well and showed both symptomatic and chemical/laboratory improvement during long-term follow-up.



## Outcome

The patient had retained stones within the cystic duct that were unable to be removed by standard balloon sweep. The SpyGlass™ DS System was used to identify a rare cystic duct variant and enabled selective cannulation. Ultimately this enabled clearance of retained stones within the cystic duct.

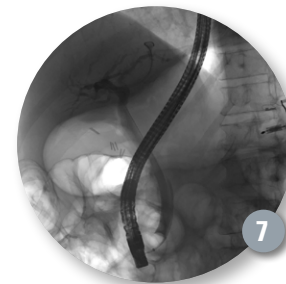

## Conclusion

The patient had retained stones within a cystic duct of unusual anatomy and low medial insertion that could not be retrieved by standard cannulation and balloon sweep. The SpyGlass DS System allowed for detection of anatomical variation in the biliary tract and selective cannulation of the cystic duct, which would have otherwise been a significant technical challenge, possibly requiring repeated procedures and surgical intervention. The SpyGlass DS System is a useful tool for delineation of biliary tract anatomy and facilitates selective cannulation, reducing procedure time and allowing for endoscopic management of disease.

**IMPORTANT INFORMATION:** These materials are intended to describe common clinical considerations and procedural steps for the use of referenced technologies but may not be appropriate for every patient or case. Decisions surrounding patient care depend on the physician's professional judgment in light of all available information for the case at hand.

Boston Scientific (BSC) does not promote or encourage the use of its devices outside their approved labeling. Case studies are not necessarily representative of clinical outcomes in all cases as individual results may vary.

All images provided courtesy of Dr. Parikh.

Results from case studies are not predictive of results in other cases. Results in other cases may vary.

All trademarks are the property of their respective owners.

**CAUTION:** Federal (USA) law restricts this device to sale by or on the order of a physician.

**CAUTION:** The law restricts these devices to sale by or on the order of a physician. Indications, contraindications, warnings and instructions for use can be found in the product labelling supplied with each device. Information for use only in countries with applicable health authority registrations. Material not intended for use in France.

Boston Scientific Corporation  
300 Boston Scientific Way  
Marlborough, MA 01752-1234  
[www.bostonscientific.com/gastro](http://www.bostonscientific.com/gastro)  
[www.EndoSuite.com](http://www.EndoSuite.com)

©2017 Boston Scientific Corporation  
or its affiliates. All rights reserved.

ENDO-493229-AA September 2017

